



Five therapies. One device.

Unified respiratory system makes care quick and easy



VS




55 lb

5 machines | Multiple circuits


18 lb

1 machine | 1 circuit

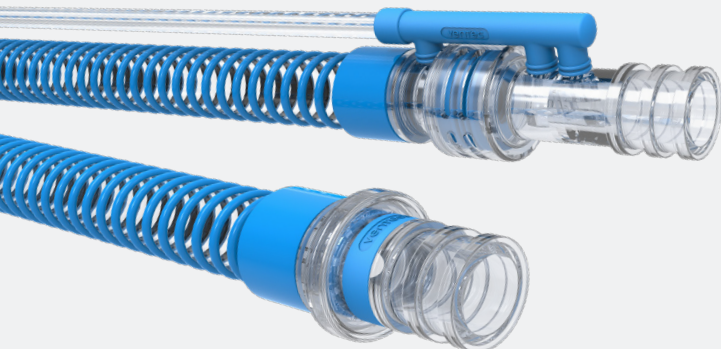

Ventec One-Circuit™

Ventec One-Circuit™ allows switching between therapies with the touch of a button. The redesigned valves make Touch Button Cough™ therapy possible and the integrated oxygen tube enables Oxygen Direct™. Rather than several cumbersome tubes connecting to separate devices, the Ventec One-Circuit™ eliminates clutter to create one easy-to-manage system for peace of mind.

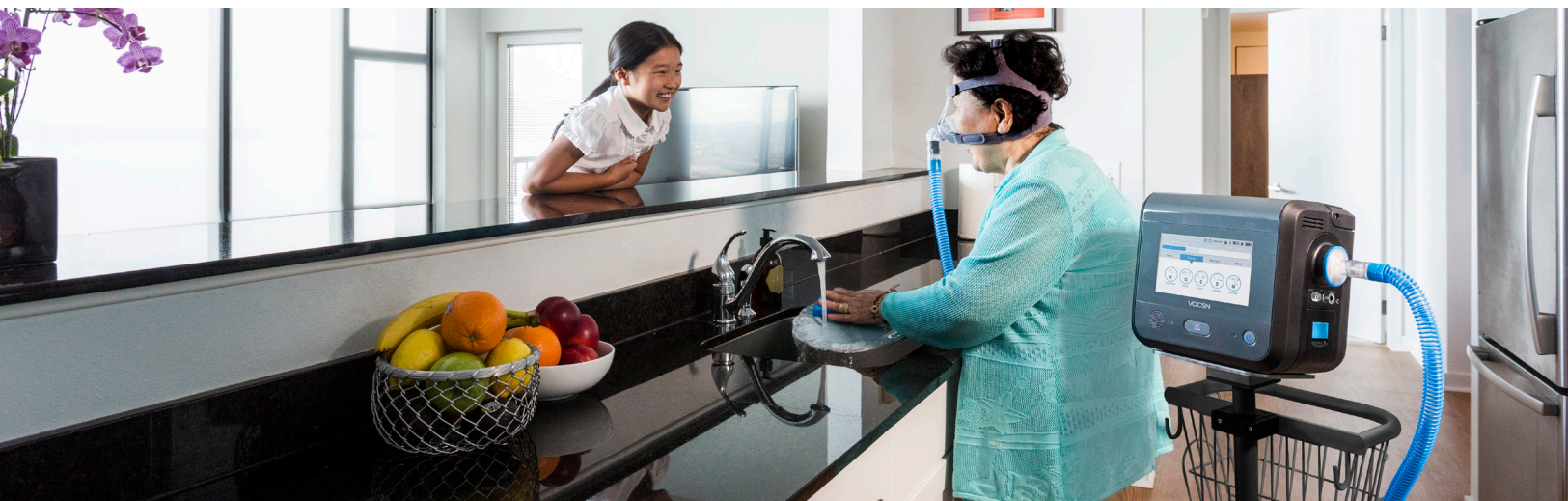

Customizable to Each Patient

Get all five therapies or just the mix of therapies needed



Ventilation

VOCSN brings next-generation technology, breathing comfort, and ease of use to ventilator patients.

Oxygen

The VOCSN internal oxygen concentrator and Oxygen Direct™ system deliver the equivalent of 6 L/min of oxygen. External oxygen sources can be connected when needed.

Cough

Touch Button Cough™ therapy can be activated in seconds, without changing the circuit.

Suction

The quiet hospital-grade suction system is easy to transport.

Nebulization

VOCSN records each medication use and turns off the nebulizer when the therapy is complete.



SIMPLE

- Ventec One-Circuit™
- Intuitive touchscreen operation
- Quick and easy setup



MOBILE

- Five therapies in one device
- 70% lighter and smaller than five devices
- Up to 9 hours of on-board battery



CARE CHANGING

- Touch Button Cough™
- Uninterrupted ventilation between therapies
- Designed to reduce risk of circuit misconnects